


waterjet[®]
FLOSSER

Oral Irrigator
Model: WaterJet Nano
WJ-360

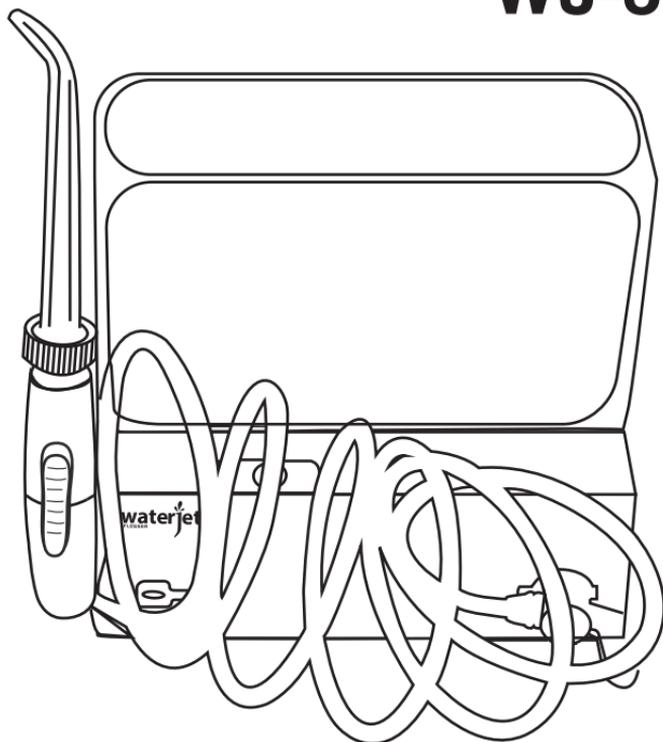


CONGRATULATION!

Thank you for selecting the WaterJet Flosser oral irrigator.

The WaterJet Flosser oral irrigator cleanses your teeth and gums, removing food particles from between your teeth and reducing bacteria that can lead to gingivitis and other periodontal diseases.

To maximize your benefit from the WaterJet Flosser oral irrigator, it is important that you read and understand the information in this instruction manual.

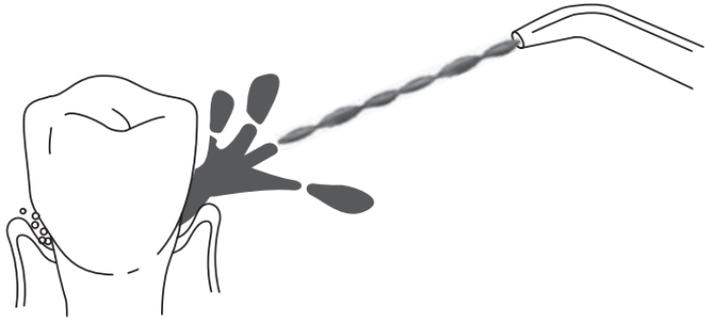

CONTENTS

Important Safeguards	4-5
Features	6
Operation	7-9
Maintenance	10
Troubleshooting	11
Warranty	12

IMPORTANT SAFEGUARDS

SAFETY:

When using electrical appliances, especially when children are present, precautions should always be followed, including, but not limited to the following:

READ ALL INSTRUCTIONS BEFORE USING.

DANGER:

To reduce the risk of electrocution:

- Always unplug this product immediately after use.
- Do not handle charger with wet hands.
- Do not place in or drop into water or other liquid.
- Periodically check the cord for damage.
- Do not use the unit while bathing.
- Do not place or store the unit where it can fall or be pulled into a tub, sink, or toilet.
- Do not reach for a unit that has fallen into water. Unplug immediately.

WARNING:

To reduce the risk of burns, electrocution, fire, or injury to persons:

- Do not plug this device into voltage system that is different from the voltage system specified on the device.
- Never operate this product if it has a damaged cord or plug, if it is not working properly, if dropped or damaged, or dropped into water.
- Use of a battery other than the special battery designed for this appliance could result in battery explosion.
- Do not direct water under the tongue, into the ear, nose or other delicate areas. This appliance is capable of producing pressures that may cause serious damage in these areas.
- Use this product only for its intended use as described in the manual or as recommended by your dental professional.
- Do not use attachments not recommended by the manufacturer
- Fill reservoir with warm water or other mouthwash or prescribed solution as recommended by your dental professional
- Never drop or insert any object into any opening or hose.

IMPORTANT SAFEGUARDS

- Keep charging unit and cord away from heated surface
- Do not use while wearing any oral jewelry. Remove prior to use.
- Do not use where oxygen or aerosol sprays are being administered.
- Do not add soda, salt or any other abrasive substance to the water in the reservoir as doing so will damage your WaterJet Flosser.
- If you have had oral or periodontal surgery in the previous twelve months, contact your dentist prior to use.
- Do not use this product for more than 8 minutes in each two-hour period.
- Do not use WaterJet Nano for more than four tanks water in each two-hour period
- Do not wash the main unit with water. This could damage the unit
- Do not dismantle the unit
- Do not submerge the unit in water
- **Closely instruct and supervise children 8 years and above and individuals with special needs in the proper use of this product**

- Children should be supervised to ensure that they do not play with the appliance
- There are no consumer serviceable parts in this appliance and it does not require maintenance


SAVE THESE INSTRUCTIONS.

FEATURES


travel case




reservoir

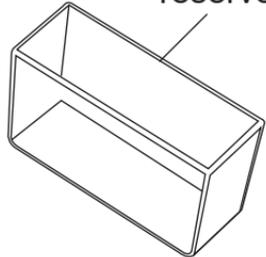


base unit

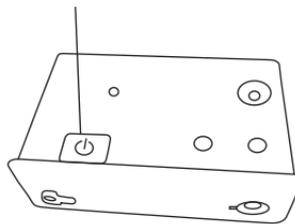

charger

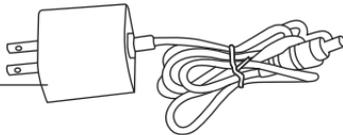

handle + hose

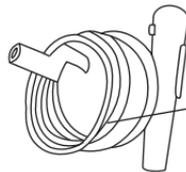

Jet Tip



Tongue Cleaner Tip



Orthodontic Tip





Periodontal Tip




OPERATION

Assembly

Insert the base of the water hose into the hole on the front of the unit. Twist the base one quarter turn to the right to lock in place

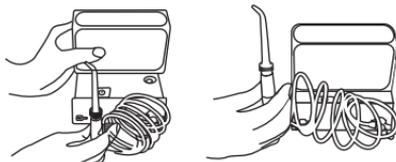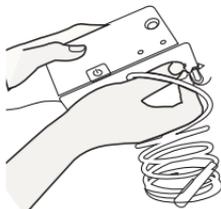


Direct the jet tip into the bottom of the sink basin and turn the on/off switch on the box of the unit to the "on" position.

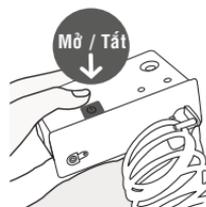

Insert jet tip firmly into the handle. When tip is correctly inserted there should be no gap between handle and the colored tip ring.

To remove tip from the handle, press the tip eject button placed on the side of the handle and pull tip from handle


Removing the tips without pressing the tip eject button can damage the handle.

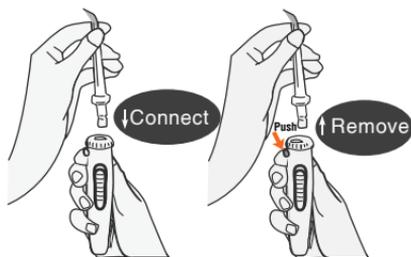

Placement in Mouth

Lean low over sink and place the jet tip in mouth. Aim the tip toward teeth.

As soon as water begins to flow from the jet tip, slide the pressure control slider to the setting that you intend to use, that is comfortable to you

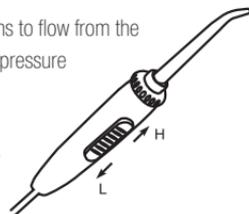

Lowest position allows you to pause completely the water flow

Direct the water flow by turning the ring at the top of the handle




Preparing the Reservoir

Remove the reservoir, and fill it with lukewarm water. Center the reservoir over the base and push down the lock in place

OPERATION

Recommended Technique

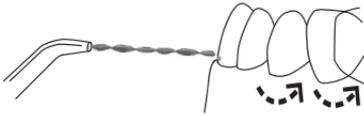

Bend over the sink and guide the jet to your teeth and gums. Start with the back teeth and work toward the front teeth. Glide tip along gumline and pause briefly between teeth. Continue until you have cleaned the inside and outside of both the upper and lower teeth. Keep your mouth slightly open to allow the water to run into the sink.

Direct the stream of water towards the gumline at 90 degrees, so that contact is made where the tooth and gums meet. Manoeuvre the flow of water from tooth to tooth remembering to pause briefly between each tooth.

You can temporarily stop the flow anytime by moving the pressure control switch at the lowest position

When Finished

Turn the unit off. **Empty any water remaining in the reservoir after use.**


If you used mouthwash or other solution, after each use, rinse the unit to prevent clogging by partially filling reservoir with warm water and running unit with tip pointed into sink until unit is empty.

First time use:


Adjust pressure control at maximum pressure, point the handle into sink and turn on the device

Charge Guide


Plug in the charger to the electrical power outlet, the LED indicator will flash green.

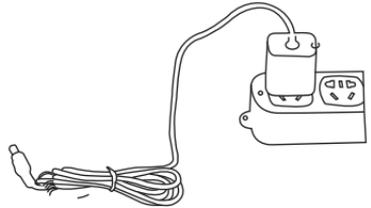

Plug the other round end of the charger into the unit's front side charging hole, the LED indicator will flash red. When red LED is on, the charging process is in progress.

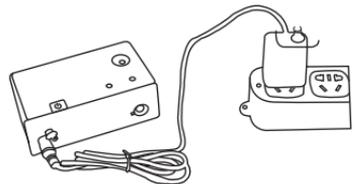

A full charge requires 15 hours. Each charge must not exceed 20 hours.

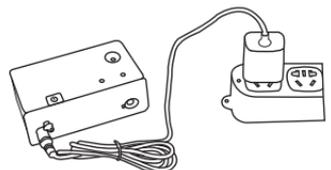

SPECIAL ACCESSORIES

Periodontal Tip


Place the periodontal tip in contact between teeth and the gumline. Lean the soft tip against a tooth at a 45-degree angle and gently place the tip under the gumline,



into the pocket. Start the irrigator by setting it at the lowest pressure and clean by gently following gumline, inserting the periodontal tip in between teeth. Continue tracing along the gumline to massage the gums and safely deliver solution into the gingival pockets below the gumline.

Orthodontic Tip

Use the Orthodontic Tip to simultaneously brush and rinse hard-to-reach areas around braces. Start on the last molar (back tooth). Place tip at 90-degree angle to

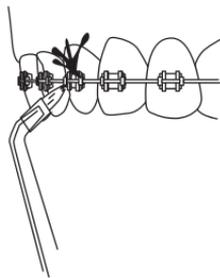

the tooth at the gumline. Gently glide tip along gumline, pausing briefly to lightly brush area between teeth and all around orthodontic bracket, before proceeding to the next tooth and outside of both the upper and lower teeth. Continue process until you have cleaned the inside.

Tongue Cleaner Tip

Lean your head over the sink. Glide the pressure button to the lowest setting. Start the irrigator and place the tongue cleaner at the back of your tongue and move it



forward to clean your tongue. Repeat to cover entire tongue. Set the pressure to the lowest setting and adjust to a higher setting for more water pressure as it feels comfortable.

MAINTENANCE

Make sure to throw out the water inside the tank, and run the appliance for 2 to 3 seconds to drain the water after usage. Otherwise bacteria could collect in the remaining water.

Do not submerge the unit in water or any other liquid.

Wipe the exterior when necessary with a soft cloth and mild non-abrasive cleanser.

Reservoir may be washed with soapy water. Before cleaning, unplug from the electrical outlet. Wash the dirt off from inside using a brush or similar, and wipe off with a cloth.

Dry the reservoir when not using for a long period.

There are no consumer serviceable parts in this appliance and it does not require maintenance

Rinse after using with mouthwash

Battery Care

To charge your WaterJet Flosser, plug the appliance into an electrical outlet. Make sure you push the switch into the "OFF" position before charging. The unit will not charge if the switch is in the "ON" position. A full charge will take at least 15 hours. Do not run the batteries down completely as this may shorten the lifespan of the batteries and will require charging the unit for 15 hours to obtain a full charge. If the WaterJet Flosser is to be stored for an extended period of time, charge

it for 15 hours before storing. Never store the WaterJet Flosser for extended periods unless it is fully charged.

Battery Disposal

This product is provided with a nickel-hydride battery, which does not contain any harmful heavy metals. At the end of its useful life, please recycle the unit at your nearest approved recycle center or disposing locations in accordance with local regulations.

TROUBLESHOOTING

PROBLEM	CAUSE	CORRECTIVE ACTION
Jet Tips and Tongue Cleaner		
Leakage between the tip and handle ring	Tip not fully engaged	Remove tip and reinstall
Periodontal Tip		
Tear in soft tip	Tip is worn	Replace with new tip
Base		
Hose to the handle leaks	Damage to the hose	Can be replaced by authorized service center or contact our company for a self service repair kit
Reservoir leaks	Missing reservoir valve or reservoir valve is upside down.	Replace reservoir valve. Turn valve right-side up
Inadequate pressure	Reservoir not positioned properly	Re-position reservoir. Press down firmly on the reservoir.
Unit doesn't start	Motor not running	Be sure outlet is functioning; test with another appliance (check electrical reset button and/or wall switch)

WARRANTY

WaterJet will repair or replace (at our option) your unit free of charge for 24 months from the date of purchase if the appliance is defective in workmanship or materials.

This guarantee does not cover: damage due to improper use, normal wear or use as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorized persons and if original WaterJet parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your purchase receipt to an authorized WaterJet distributor/customer service centre in your country. Purchaser is responsible for all postage and freight charges and any other expenses incidental to claiming against this Warranty.



WaterJet

600 W Santa Ana Blvd Ste 114A PMB 224

Santa Ana, CA 92701, USA

Made in China

www.waterjet-flosser.com