


waterjet®
FLOSSER PRO
Oral Irrigator

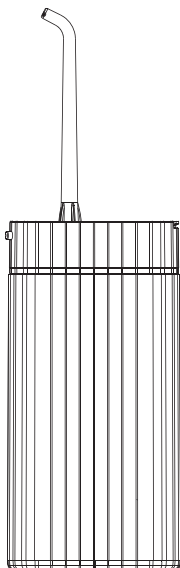

Model:

Cordless Mini Blue
WFP-440


Cordless Mini White
WFP-400

FC CE RoHS IPX7

BEFORE YOU BEGIN

The WaterJet Flosser Pro - Cordless Mini oral irrigator helps cleanse your teeth and gums, removing food particles from between your teeth and reducing bacteria that can lead to gingivitis and other periodontal diseases.

To maximize the benefits of the WaterJet Flosser Pro - Cordless Mini oral irrigator, it is important that you carefully read and understand the information provided in this instruction manual.

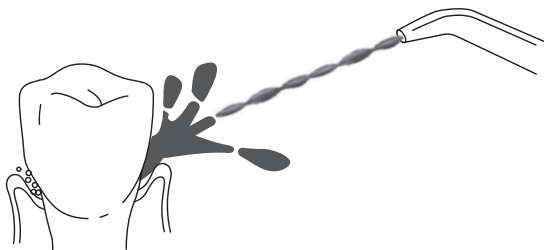

FRIENDLY REMINDER

As your gums are not yet accustomed to the stimulation, mild bleeding may occur during initial use. Please start with the lowest mode first. This bleeding should stop within a week. However, if the bleeding persists for more than two weeks, it may be a sign of a gum condition (e.g., periodontitis). In that case, please stop using the product and consult your dentist.

Safety Instructions	4
Parts Overview	5
Instructions for Use	6 - 7
Mode Indication	8
Technical Parameters	8
Product Maintenance	9
Fault Diagnosis	10
Warranty	11

SAFETY INSTRUCTIONS

DANGER

If you do not know the operating process and safety instructions, please do not rush to use the product to avoid hurts caused by mis-operation.

Do not disassemble the product or change the parts of the product to avoid damage on the product in case may caused in-valid warranty.

Do not use above 50°C hot water into the water tank for cleanning mouth to avoid scald.

If the power cable or plug is damaged, please replace it before using.

Do not put the product into the fire or the heating, and do not charge, use or leave it at the high temperature.

WARNING

Do not use this product for any other purpose other than oral cleaning.

Do not push the jet tip hard on the gums or teeth to avoid hurting on the gums.

Do not immerse this product into water for a long time to avoid incorrect operation or cause fire, electric shock and other dangers.

Do not damage, destroy or modify the power cable, and do not stretch, bend or twist the power cable. Do not press heavy objects on the power cable or get it caught between objects, and keep the product and cable away from hot objects.

Do not store the product into a hot or humid places, and avoid direct sunlight.

Do not charge the product in hot or humid places, such as the bathroom or toilet.

If not in use of the product for a long time, please unplug the power plug and charging plug to avoid fire or danger.

NOTICE

There are no replacements or repairs for this product. If repair required, please return it to the seller.


Do not use without a nozzle.

This product is not suitable for children and infants, unless they are supervised or guided by their safety responsible supervisor on how to use this product. Keep out of reach of children and infants. Do not use it as a toy.


If you have periodontal disease or have had oral surgery in the past 2 months, please consult your dentist before using.

Do not add mouthwash into the water tank which will cause corrosion and the tank will be adhered to the device body can result in tank opening failed.

PARTS OVERVIEW

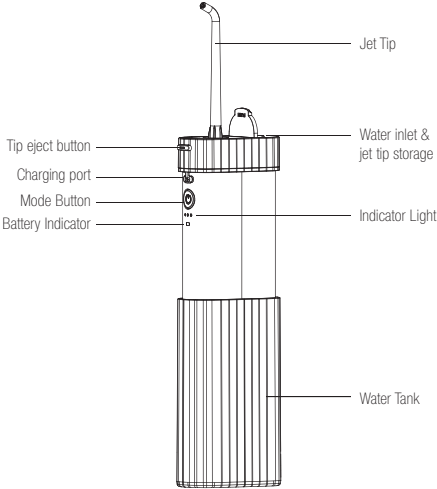


Accessories:

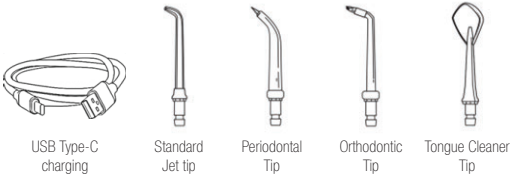

INSTRUCTIONS FOR USE


1. Attach the jet tip to the unit

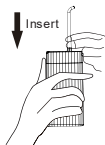


2. Pull the water tank to the very bottom to avoid water leakage and then open the tank cover to fill in with water



3. Put the Jet tip into the mouth and press the power button to start flossing

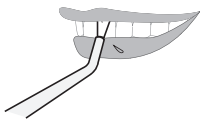


4. Long press the power switch button for 2s to switch the mode, short press to turn on/off

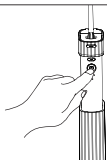

5. Storage after using: Remove the jet tip by pressing the tip eject button, open the water tank cover and store the tip.




6. Regular cleaning: Unattach the tank box and clean it

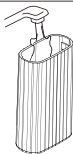

INSTRUCTIONS FOR USE

TIP USAGE

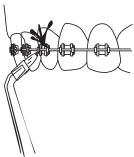 A line drawing showing a side view of a person's head bent over a sink. A water flosser tip is directed at the teeth and gums. Arrows indicate the water stream hitting the gumline.	Standard Jet Tip Bend over the sink and aim the jet tip at your teeth and gums. Press the power button to start the water stream. Keep your mouth slightly open to allow the water to flow into the sink. Direct the water stream towards the gumline at a 90° angle, ensuring contact where the teeth and gums meet. Move the tip from tooth to tooth, pausing briefly between each one.
 A line drawing showing a top-down view of a tooth. A water flosser tip is inserted between the tooth and the gumline at a 45-degree angle. Arrows show the water stream flowing into the periodontal pocket.	Periodontal Tip Place the periodontal tip between the teeth and along the gumline. Lean the soft tip against a tooth at a 45° angle and gently insert it under the gumline, into the pocket. Start the irrigator at the lowest pressure setting. Gently follow the gumline, inserting the tip between each tooth. Continue tracing along the gumline to massage the gums and deliver the water safely into the periodontal pockets below the gumline.
 A line drawing showing a side view of a person's teeth with braces. A water flosser tip is directed at the gumline between the teeth and around the brackets. Arrows indicate the water stream cleaning these areas.	Orthodontic Tip Use the Orthodontic Tip to brush and rinse hard-to-reach areas around braces, crowns, bridges, and implants. Start at the last molar (back tooth). Place the tip at a 90° angle to the tooth at the gumline. Gently glide the tip along the gumline, pausing briefly to clean between the teeth and around orthodontic brackets before moving on to the next tooth. Repeat for both the upper and lower arches, cleaning the inside and outside of each tooth.
 A 3D illustration of a person's mouth with the tongue out. A water flosser tip is shown cleaning the surface of the tongue. Small bubbles indicate the water stream.	Tongue Cleaner Tip Lean your head over the sink. Start the irrigator at the lowest setting. Place the tongue cleaner at the back of your tongue and move it forward to clean the surface. Repeat as needed to cover the entire tongue. You may gradually increase the pressure setting for a deeper clean, as long as it feels comfortable.


MODE INDICATION

Mode	Mode switch	Mode indication color	Target user	Water pressure
Gentle Mode	Press the switch to start flossing	Green	Gentle water pressure. Gentle cleaning. Suitable for sensitive teeth and first-time users.	60-80 PSI
Pulse Mode	Press the switch and hold for 2s to switch to pulse mode	Blue	Variable frequency pulse pressure. Alternates between strong & weak modes to deeply clean and massage the gums.	60-110 PSI
Strong Mode	Press the switch and hold for 2s to switch to strong mode	White	Strong water pressure for efficient cleaning, suitable for daily use by most people.	120-140 PSI
Battery indicator	1.The indicator turns red when run out of power. 2.The indicator flashes red when charging. 3.When fully charged, the indicator turns green.			

TECHNICAL PARAMETERS

Product name	WaterJet Flosser Pro Cordless Mini	Product model	WFP-400 (White) WFP-440 (Blue)
Jet tip diameter	0.6mm	Pulse frequency	1600 times/minute
Water pressure	60-140PSI	Water tank capacity	140ml
Battery capacity	1100mAh	Noise	≤72dB
Input power	5V=1A	Rated power	8W
Charging port	Type-C	Charging time	About 4h
Waterproof level	IPX7	Water pressure mode	3 modes
Product size (Height*width*thickness)	127*33*63mm	Net weight	198g ± 5g

PRODUCT MAINTENANCE

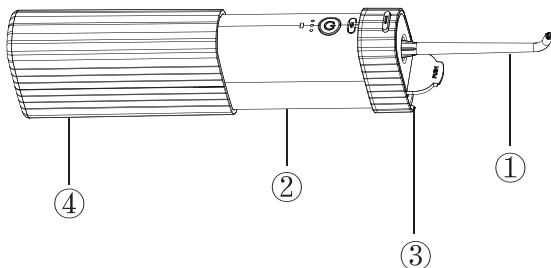

1. Nozzle (tip)

- Press the tip eject button and remove the nozzle (tip), wash it with water and then wipe dry with a cloth.
- Dentists recommend changing the nozzle (tip) every 6 months.
- If the top of the nozzle (tip) is deformed, replace it with a new one to avoid malfunction.

2. Product body

- Wipe off any dirt using a soft cloth dampened with liquid soap, then dry with a clean cloth. If the product body is dirty, wash it with water.
- Do not use thinner, gasoline, alcohol, or other harsh solvents to clean the product. Do not wash with hot water above 50°C. After cleaning, please dry the device thoroughly.

3. Water inlet

- Wipe off the dirt, otherwise it will cause a leakage

4. Water tank

- Open the tank and clean it regularly.
- Clean the inside using a brush or similar tool, then wipe dry with a clean cloth.
- If not using the product for a long time, please keep the water tank dry.

FAULT DIAGNOSIS

Problem	Possible reasons	Solution
No working	Newly purchased or not used for more than three months	Please charge it for more than 4 hours
Low water pressure when working	Something inside the Jet tip and blocked the Jet tip into deformation	Replace the Jet tip
Battery only can be used for a short time when fully charged	The battery is reaching the end of its usable life	Contact your seller
	Insufficient charging times	Recharge the battery for at least 4 hours before using
The charging indicator light does not light up when charging	The power plug is not properly inserted	Please connect the power correctly
	Product charging circuit is faulty	Please contact your seller

WARRANTY

WaterJet will repair or replace (at our option) your unit free of charge for 12 months from the date of purchase if the appliance is defective in workmanship or materials.

This warranty does not cover: damage due to improper use, normal wear or use as well as defects that have a negligible effect on the value or operation of the appliance. The warranty becomes void if repairs are undertaken by unauthorized persons and if original WaterJet parts are not used.

To obtain service within the warranty period, hand in or send the complete appliance with your purchase receipt to an authorized WaterJet distributor/ customer service center in your country. Purchaser is responsible for all postage and freight charges and any other expenses incidental to claiming against this Warranty.



WaterJet

600 W Santa Ana Blvd Ste 114A PMB 224
Santa Ana, CA 92701, USA

www.waterjet-flosser.com



Dongguan Mingchuang Youshang
Electronic Technology Co., Ltd.

Bldg. 5, Haiyong Tech Park, Fenggang,
Dongguan, Guangdong, China

Made in China