


waterjet®
FLOSSER

Oral Irrigator

**Model:
Cordless Maxi
WJ-900**

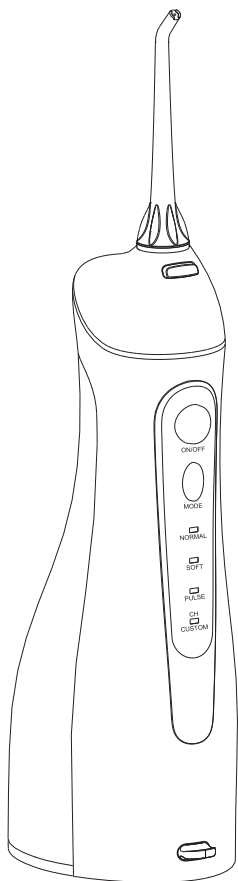


CONGRATULATION!

Thank you for selecting the WaterJet Flosser oral irrigator.

The WaterJet Flosser oral irrigator cleanses your teeth and gums, removing food particles from between your teeth and reducing bacteria that can lead to gingivitis and other periodontal diseases.

To maximize your benefit from the WaterJet Flosser oral irrigator, it is important that you read and understand the information in this instruction manual.

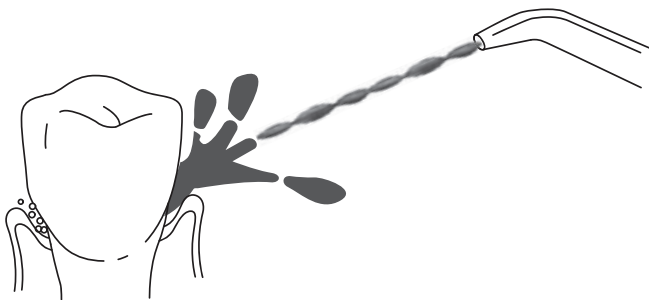

CONTENTS

Important Safeguards	4-5
Operation	6-8
Maintenance	9
Disposal	10
Troubleshooting	11
Warranty	12

IMPORTANT SAFEGUARDS

All information in this section is provided to instruct you in the correct and safe operation of this device so that you can avoid injuries to yourself and others as well as damage to property.

SAFETY:

When using electrical appliances, especially when children are present, precautions should always be followed, including, but not limited to the following:

READ ALL INSTRUCTIONS BEFORE USING.

DANGER:

To reduce the risk of electrocution:

- Toddlers, infants, people who are unable to operate should not allowed use this appliance.
- Keep out of children's reach
- Never alter the appliance. Also, do not try to disassemble or repair it.
- Do not handle charger with wet hands.
- Never alter the appliance or try to disassemble, repair it.
- Do not use the unit while bathing.
- Do not place or store the unit where it can fall or be pulled into a tub, sink, or toilet.
- Do not reach for a unit that has fallen into water. Unplug immediately.

WARNING:

- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children should be supervised to ensure that they do not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
- The appliance is only to be used with the power supply unit provided with the appliance. Use of other power supply can cause explosion
- The battery must be removed from the appliance before it is scrapped.
- The appliance must be disconnected from the supply mains when removing the battery
- The battery is to be disposed of safely
- Regarding how to remove the battery, please to see "Disposal" in page 8.
- Do not fill the tank with water warmer than 40 degree Celsius
- If the external flexible cable or cord of this transformer is damaged, it shall be replaced by the manufacturer or his service agent to avoid a Hazard.
- Keep charging unit and cord away from heated surface

IMPORTANT SAFEGUARDS


SAVE THESE INSTRUCTIONS.

- Do not use attachments not recommended by the manufacturer
- Do not use if you lost the nozzle
- Do not submerge the unit in water
- Do not use for any purpose other than oral cleaning
- Keep nozzle away from children. The nozzle is small enough to be swallowed and get lodged in the throat
- Do not fill mouthwash in the water reservoir to avoid cracking
- Certain types of mouthwash may damage the equipment by causing the housing or water reservoir to crack. Therefore, do not use mouthwash in the water reservoir.
- Do not add soda, salt or any other abrasive substance to the water in the reservoir as doing so causes damage.
- Do not wash the main unit with water. This could damage the unit
- People suffering from tooth or gum pain may not use this device and should consult their dentist
- Do not press the nozzle too firmly against your teeth and gums to avoid injuries.

OPERATION

Assembly


Attach nozzle onto the unit

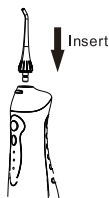

Open the water tank cap to fill the tank with water which temperature not exceed 40°C, then close the cap



Or slide the water tank down off the unit, fill the tank with water, then attach the tank onto the unit

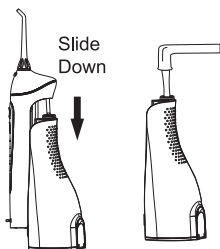

Press the "MODE" key to select the mode you desire. During irrigation, the working mode can be changed. Cordless Maxi provides the largest range of pressures of all Cordless Oral Irrigators; from strong pressure for efficient cleaning to very low pressure for sensitive and inflamed gums.

NORMAL: Pressure 90 PSI - for regular efficient cleaning.


SOFT: Pressure 40 PSI - for gentle cleaning of sensitive gums.

PULSE: Pressure 40 PSI with regular pulsation - for cleaning and massaging the gums.

CUSTOM: Pressure between 10 PSI to 90 PSI with 8 levels- Adjust precisely to your ideal pressure level (New units from factory have setting at 10 PSI, ideal for patient with very sensitive gums naturally or after a dental procedure)

Custom setting procedure:

- 1) Press the "ON/OFF" key to turn unit ON (without water)
- 2) Repeatedly press the "MODE" key until the light arrive on "CUSTOM"
- 3) Once on "CUSTOM" press and hold the "MODE" key for 3-5 seconds until the blue lighting above "CUSTOM" begins to blink.
- 4) After the light begins blinking, let go of the "MODE" key and immediately press it again repeatedly until the desired pressure is reached. At each pressure you should hear the sound of the motor change between 8 levels of intensity from lowest pressure (low noise) to highest pressure (high noise). The "CUSTOM" light should always be blinking throughout this step.
- 5) At the level of pressure you want to keep saved, stop pressing the "MODE" key until the "CUSTOM" light stops blinking and stays on continuously - Your pressure level for "CUSTOM" mode is now saved.

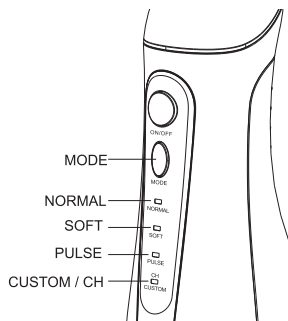


OPERATION

Hold the unit in an upright position with nozzle tip point toward your teeth.

Keep your mouth opened slightly so that the water can flow out.


Press "On/Off" key to start irrigating.

Press "On/Off" key and then press the nozzle release button to replace other nozzle when you need.

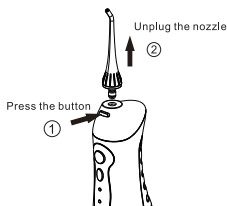

Once irrigating complete, press "On/Off" key again to turn off the device

Recommended Technique

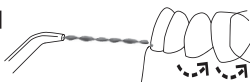


Directing the water jet flow at a 90-degree angle to the teeth and gums

Directing the water jet flow slowly along the teeth

Directing the water jet along the gum line with the tip at a 90-degree angle to the gum tissue

Cleaning between the teeth and around any braces, crowns, or bridgework thoroughly

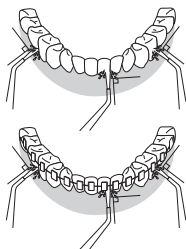


When Finished

Turn the unit off. Open the water tank cap and dispose of the water in the tank.

Turn on the unit and drain the water inside the appliance


Wipe the device with dry cloth ASAP before store it.

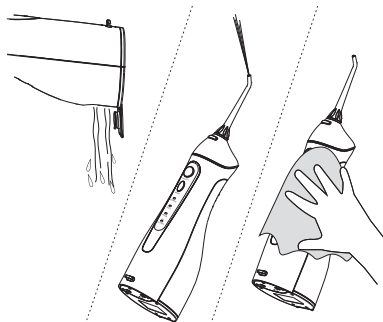

OPERATION

Charge Guide

Take the charger base cover down and insert the appliance plug securely into the appliance socket (make sure no liquid remaining inside the socket before insert the appliance plug).

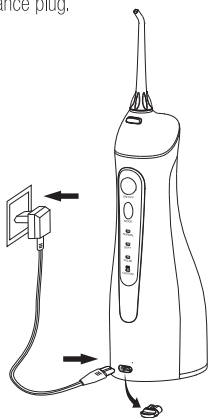

Insert the charging adaptor securely into the power outlet

The "CH" indicator blinks when the device is charging, and it will be steadily on after the device is full charged

Recover the cover onto charger base again after charging finished

Note:

If the "CH" indicator is blinking during or after an irrigating, that means the battery need to be recharged. A full recharging will take about 4 hours.

MAINTENANCE

Use only water or neutral cleaning agent to clean the device. Never use caustic or abrasive agents, since they may damage the device.

Remove the dirt by wiping using a cloth with liquid soap and wipe off using a dry cloth

Do not submerge the unit in water or any other liquid.

Do not wash using hot water of over 50 degree Celsius

Keep the device away from high temperature area and avoid direct sunlight

Wash nozzle with water and wipe off with soft cloth.
Replace the nozzle every 6 months

Do not bend, pull or twist the water suction tube

Insert a screwdriver into the slot located at the bottom of the handle and turn counter-clockwise to release the bottom cap




ISO9001 & ISO13485 Certified Manufacturer




Printed in China

DISPOSAL


The rechargeable battery inside your appliance is easy to remove for recycling. At the end of the product's life and prior to disposal, please remove the rechargeable battery from the handle by following the instruction below. During below operation, be sure to protect your eyes, hands, fingers, and the surface on which you work.

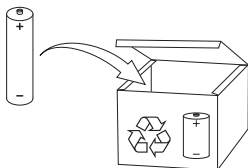


Remove 2 pcs rubber feet on the bottom the handle. Insert a screwdriver into the slot located at the bottom of the handle and turn counter-clockwise to release the bottom cap.

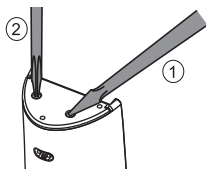


Slide the water tank down off the device and unscrew the two screws beside of the suction tube.

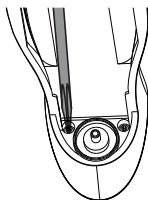

Open the top cover and put the inner out of the unit, then unscrew all the screws of the inner

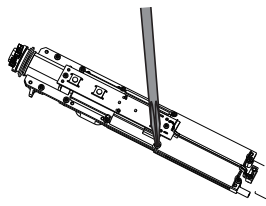

Cut off the battery line and pull out the batteries from the inner

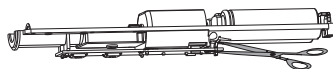

If you have recycling questions, please contact your local waste management office

TROUBLESHOOTING

PROBLEM	CAUSE	CORRECTIVE ACTION
It does not work	You've just purchased the product or it has been left unused for more than 3 months	Charge the unit for at least 8 hours
It only works for a few minutes even after charging	The battery has reached the end of its lift cycle	Contact service center support
	Charging time was insufficient	Charge the unit for at least 8 hours
Water pressure is very low	There is dirt stuck in the nozzle	Replace the nozzle
	The nozzle is deformed	
	There is dirt stuck in filter	Clean the filter
Water does not come out	The tank is empty	Fill up tank
	The appliance is tilted too much while in use	Use with the appliance in a vertical state

WARRANTY

WaterJet will repair or replace (at our option) your unit free of charge for 12 months from the date of purchase if the appliance is defective in workmanship or materials.

This guarantee does not cover: damage due to improper use, normal wear or use as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorized persons and if original WaterJet parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your purchase receipt to an authorized WaterJet distributor/customer service centre in your country. Purchaser is responsible for all postage and freight charges and any other expenses incidental to claiming against this Warranty.



WaterJet

600 W Santa Ana Blvd Ste 114A PMB 224

Santa Ana, CA 92701, USA

Made in China

www.waterjet-flosser.com