


waterjet®
FLOSSER PRO
Oral Irrigator


Model:
Slim
WPS-300

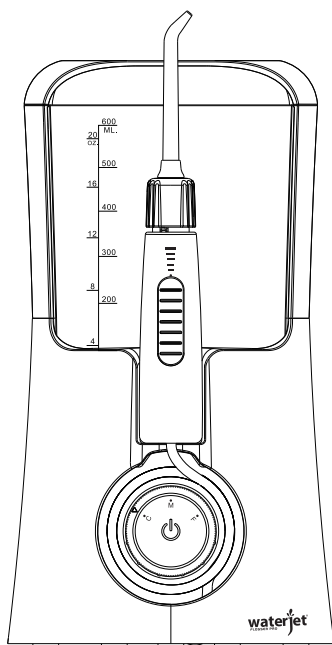

CONGRATULATION!

Thank you for selecting the WaterJet Flosser oral irrigator.

The WaterJet Flosser oral irrigator cleanses your teeth and gums, removing food particles from between your teeth and reducing bacteria that can lead to gingivitis and other periodontal diseases.

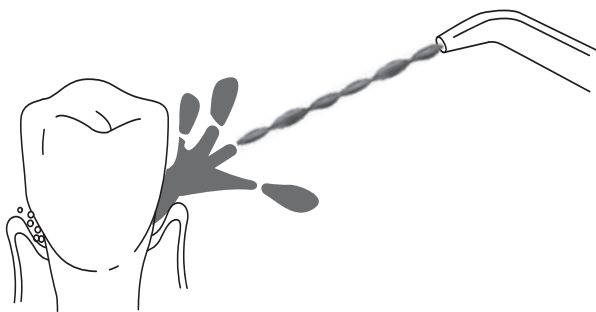

To maximize your benefit from the WaterJet Flosser oral irrigator, it is important that you read and understand the information in this instruction manual.

CONTENTS

Important safeguards	5-6
Features	7
Operation	8-12
Maintenance	13
Troubleshooting	14
Disposal	15
Warranty	16

SAFETY:

Please read the safety instructions before using this product. The safety instructions provided here help you to use this product correctly, and also to prevent you from injuring yourself or others

When using electrical products, especially when children are present, the basic safety precautions should always be followed.

READ ALL INSTRUCTIONS BEFORE USING.

DANGER:

To reduce the risk of electrocution:

- Toddlers, infants, children should not be allowed to use this device
- Please keep it out of reach of children
- Do not use if the power cord or plug was damaged
- Always unplug the device after using
- Do not touch plug with wet hand
- Do not place the device in water or other liquids
- Do not place the device where it may fall or be pulled into a tub or sink
- Do not use the device while bathing
- Do not touch the product which has fallen into water. Unplug immediately

WARNING:

To reduce the risk of burns, electrocution, fire, or injury to persons:




- Cleaning and user maintenance shall not be made by children without supervision
- This appliance can be used by children aged from 8 years and above and people with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved
- Children should not play with the appliance
- The appliance is only to be used with the power supply unit provided with the appliance
- Keep the device away from heated surface
- Do not use other tips except recommended parts
- Please do not use the device if you lose the nozzle
- Keep nozzle away from small children. The nozzle is small enough to be swallowed and can become lodged in the throat

IMPORTANT SAFEGUARDS



- Do not direct water under the tongue, into the ear, nose or other delicate areas
- Do not use for any other purpose except oral cleaning.
- Do not fill the tank with warm water over than 40 degree C.
- Certain types of mouthwash may damage the device by causing the water tank to crack. Therefore, do not use mouthwash in the water tank
- People who experience tooth or gum pain may not be able to use this device and should consult their dentist
- Do not press the nozzle too close against your teeth or gums, doing so may injure your gums
- For household use only, not for outside use
- Please fill the tank with pure water

SAVE THESE INSTRUCTIONS.

FEATURES

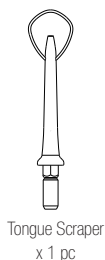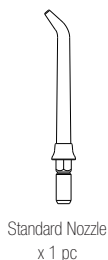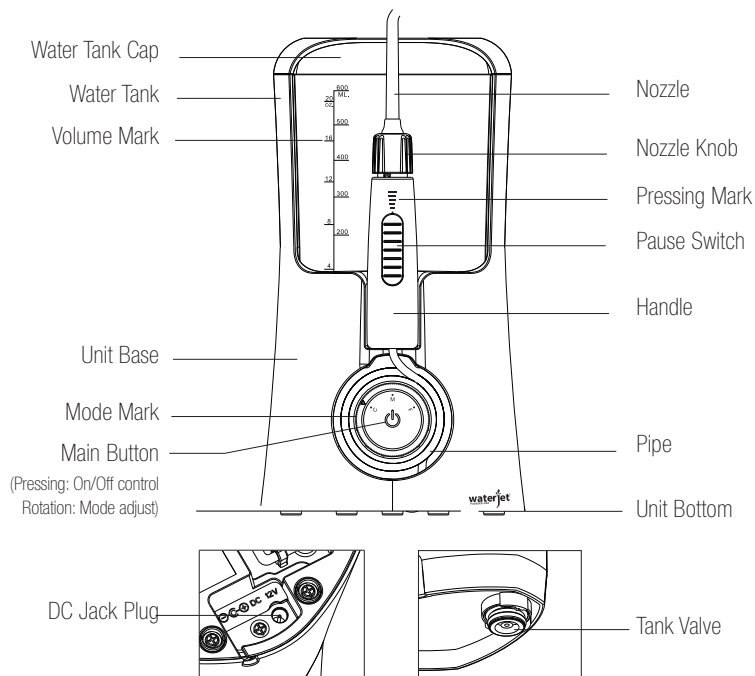

Replace every 6 months

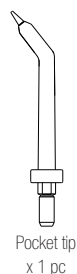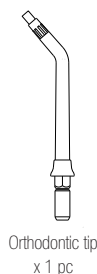

Replace every 3 months

Getting started

Place the device on a flat surface without any sharp objects and make sure the cord is dry and has no damage. The device is powered by a special adaptor with DC jack connector.

(Input: 100-240V~, 50/60Hz; Output: 12V, 1.5A)

Insert the plug into the DC Jack on the bottom of base and plug the adaptor cord into an electrical outlet (Fig.1)

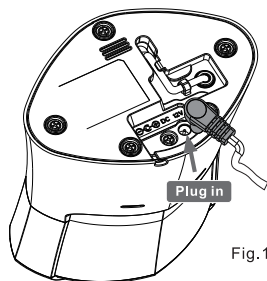


Fig.1

Preparing the Reservoir

Take the water tank cap off

Take the water tank off to fill with proper liquid, no warmer than 40 degree C, then replace the tank onto the top of main unit and press down firmly (Fig.2-1 & Fig.2-2)

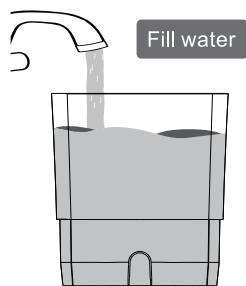


Fig.2-1


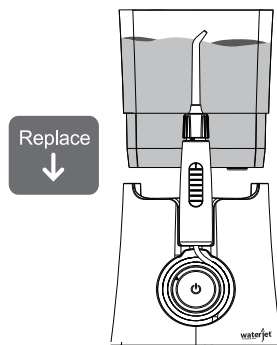


Fig.2-2

Inserting and Removing Nozzle


Inserting Nozzle

Press down the knob and insert tip into the center of the handle, also try to pull the tip out gently to confirm the tip is locked firmly. (Fig.3)

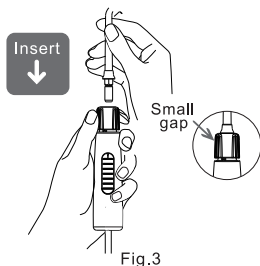

Removing Nozzle


Press the Release Button on the top end of the handle, then pull the nozzle out lightly by hand (Fig.4)

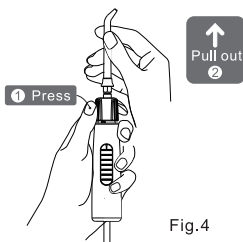



How to use the WaterJet Flosser

Do not power on the device before placing the nozzle in your mouth. First use, we suggest you choose the lowest pressure, gradually increase the water pressure over time to the setting you need.

- 1) Hold the handle and point the nozzle to sink, then press down the Main Button to power the device on (Fig.6)

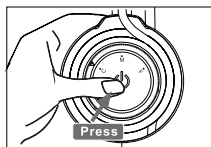

- 2) Rotate the Main Button to change mode, after setting, push the handle switch to adjust water pressure (Fig.7-1 & Fig.7-2)

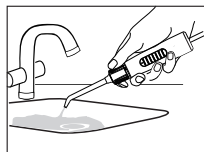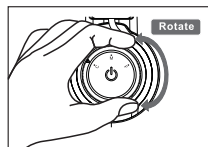

OPERATION

3) Adjust the direction of nozzle by rotating the nozzle knob

4) Bend over to the sink and place the tip in your mouth, and then push the Pause Switch to "1" position to enjoy irrigating (Slightly close your mouth to avoid the water splashing). (Fig.8)



Fig.8

5) Position the WaterJet Flosser flow at 90-degree angle to the teeth and gums. (Fig.9)

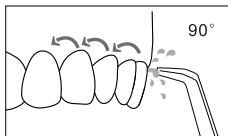


Fig.9

6) It has 3 different optional modes: Clean mode, Massage mode and Fresh breath mode

7) You want to use Clean mode, rotate switch to "C". In Clean mode water flow will short pause at 30 seconds, which reminds you to clean the other arch of your mouth and another short pause at 1-minute to remind you the recommended usage time is up. The unit will continue to work till 3-minute cycle is complete or you turn it off. (Fig.10)

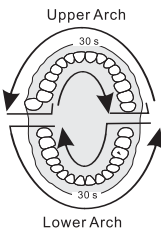


Fig.10

8) If you rotate switch to "M", Massage mode can give your teeth good massage experience. The unit will continue to work till cycle is complete or you turn it off. (Fig.11)

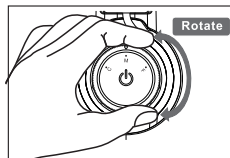


Fig.11

9) If you want to use Fresh breath mode, rotate switch to "F". In this mode, you can use 0.5% salt water or 8% mouthwash (50ml mouthwash and other water), it can give you fresher experience.

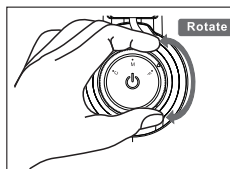


Fig.12

The unit will continue to work till cycle is complete or you turn it off. After using, you can use clean water to wash it (Fig.12)

How to end use

You can end irrigating through the below methods:

1. Auto shut-off turns off after 3-minute cycle is complete
2. Turn it off by pressing the Main Button
3. Unplug the adaptor directly when in any emergency.

When Finished

Take off the water tank and dispose of water in it

Push the slide switch to the highest pressure setting, and power on to drain the water inside the device

Turn off the device and wipe it with a soft dry cloth (Fig.13)

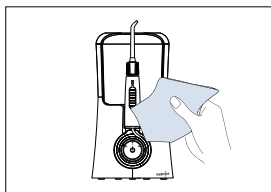


Fig.13

Press down the Release Button, and remove the nozzle, then place it into the nozzle case (Fig.14)

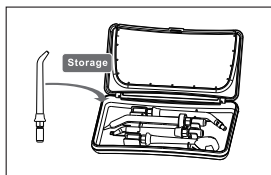


Fig.14

Unplug the adapter.

SPECIAL ACCESSORIES

Tongue Scraper

When you want to clean your tongue, please gently place the tongue scraper onto the center midline of your tongue surface about halfway back with light water pressure setting first, and gradually increase the pressure as needed. (Fig.15)

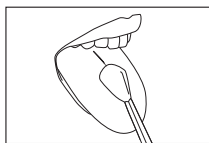


Fig 15

Orthodontic Tip

The orthodontic tip is specially designed for people wearing braces: It will clean the hard to reach areas around and in braces

Aim water stream at 90 degree to the tooth at the gum-line, start on the last molar (back tooth)

Slowly slide the tip along with the gum-line, and pausing briefly to use the brush bristles to clean the area between teeth and all around orthodontic braces

Repeat until you have finished cleaning the inside and outside of all your teeth (Fig.16)

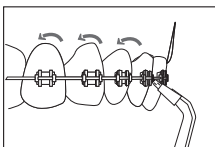


Fig.16

Pocket Tip

The pocket tip is specially designed for the cleaning of periodontal pockets. To use the pocket tip, it is a **MUST** to set the device to the lowest pressure setting

Place the soft tip against a tooth at 45-degree angle and gently place the tip under the gum-line, then insert into the pocket (Fig.17)

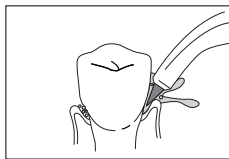


Fig.17

Turn the unit on and continue tracing along the gum line until all periodontal pockets are treated (Fig.18).

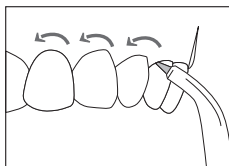


Fig.18

Note:

This tip must be used only with the lowest water pressure to extend life of soft tip as it is easily broken

Use only water or neutral cleaning agent to clean the device. Never use caustic or abrasive agents (e.g., vinegar or lime scale removers, since they may damage the device. Do not wash using hot water of over 40 degree C. Keep the device away from high temperature area and avoid direct sunlight.

Main Unit

Remove the water tank and drain the remained water inside the device.

Clean by using a soft cloth with liquid soap and dry with a soft dry cloth.

Do not submerge the main unit in water to clean

Handle

Wash with water and dry up with a soft cloth

Do not bend, pull or twist the water hose

Water Tank

Wash with water and wipe with a soft dry cloth

Make sure the water tank is inverted and firmly fixed to store

Nozzle

Wash with water and wipe with a soft dry cloth

Place the nozzle in the storage compartment

Replace once every 6 months for Standard nozzle and Tongue scraper, and once every 3 months for other tips

Charging Adaptor

Wipe the metal part of the power plug with a dry cloth once every 6 months

TROUBLE SHOOTING

PROBLEM	CAUSE	CORRECTIVE ACTION
Unit does not work	Motor not running	Make sure the adapter is firmly connected to the outlet and the unit
		Make sure the outlet and unit power is on
Water Pressure is very slow	Dirt stuck in the nozzle	Clean or replace with a new nozzle
	The nozzle is deformed	Replace with a new nozzle
	Dirt stuck in the tank valve	Clean the valve and fill the tank with pure water
Water does not come out	The tank is empty	Fill the tank up
	Water pressure under the lowest setting	Slide the switch to adjust pressure
	Tank valve not open	Press down the tank, or reassemble the tank valve
Tank leaks	The valve was lost	Replace with a new tank valve
	The valve is not installed in place	Take off the valve and reinstall again
Tear in Pocket tip head	Tip head is worn	Replace with a new nozzle

DISPOSAL

Do not disassemble the device at any time other than when disposing of it

Doing so earlier can destroy the product

The appliance must be disconnected from the power supply when disassembling the device



015

ISO9001 & ISO13485 Certificated Manufacturer



WARRANTY

WaterJet will repair or replace (at our option) your unit free of charge for 24 months from the date of purchase if the appliance is defective in workmanship or materials.

This guarantee does not cover: damage due to improper use, normal wear or use as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorized persons and if original WaterJet parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your purchase receipt to an authorized WaterJet distributor/customer service centre in your country. Purchaser is responsible for all postage and freight charges and any other expenses incidental to claiming against this Warranty.



WaterJet

600 W Santa Ana Blvd Ste 114A PMB 224
Santa Ana, CA 92701, USA

www.waterjet-flosser.com



Shenzhen Risun Technology Co., Ltd.
Building A, No. 6 of Xinmu Road,
Xinmu Community, Pinghu, Shenzhen,
China

Made in China