




**Rechargeable
Toothbrush
Model:
Smart Sonic
SW-800**

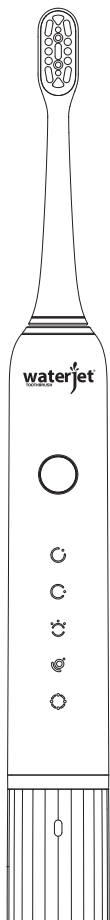


CONGRATULATION & CONFIGURATION

Thank you for selecting the WaterJet Smart Sonic Toothbrush.

The WaterJet Toothbrush cleanses your teeth and gums, removing food particles from between your teeth and reducing bacteria that can lead to gingivitis and other periodontal diseases.

To maximize your benefit from the WaterJet Toothbrush, it is important that you read and understand the information in this instruction manual.

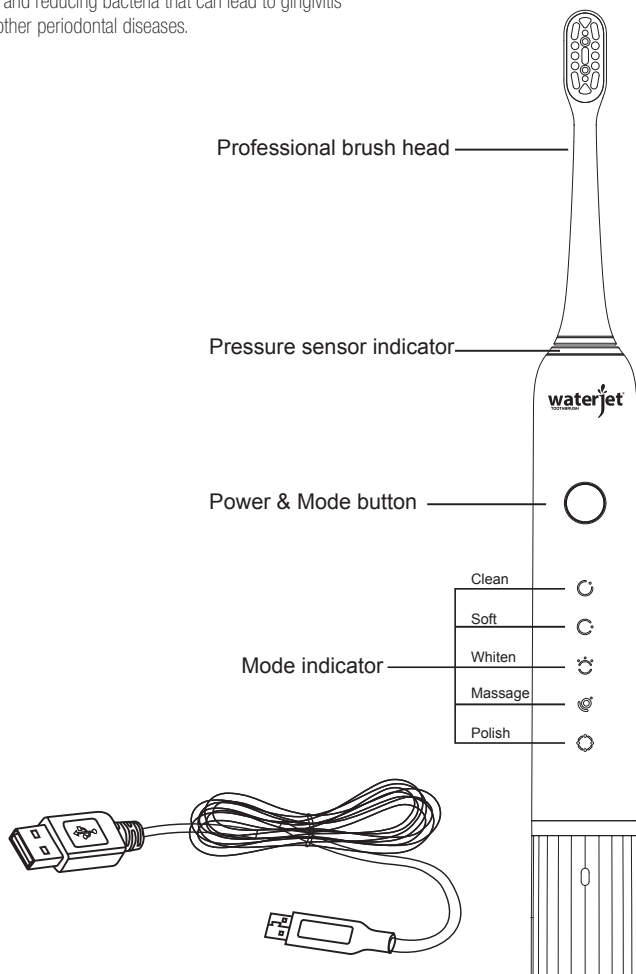

CONTENT

Important Safeguards	4-5
Getting started	6
How to use	7-8
Maintenance	9
Disposal	10
Troubleshooting	11
Warranty	12

IMPORTANT SAFEGUARDS

The Smart Sonic toothbrush is safe and effective when used according to directions. Like any electrical appliance, however, basic safety precautions should be followed, especially when children are present.

READ ALL INSTRUCTIONS BEFORE USING.

DANGERS-WATER

To reduce the risk of burns, electrocution, fire or injury to persons:

1. Do not place or store the charger base where it can fall or be pulled into a bathtub, sink or toilet.
2. Do not place or drop the charger into water or other liquid.
3. Do not reach for the product that has fallen into water or other liquid. Unplug immediately.
4. Do not use while bathing.
5. Unplug the charger base before cleaning and assure that it is dry before plugging it into an electrical outlet.

DANGERS-VOLTAGE

To reduce the risk of burns, electrocution, fire or injury to persons:

1. Do not plug the charger base into an outlet with a voltage other than specified on the bottom of the charger.
2. This product is designed to operate within a range of 100-240 volts.
3. Voltage converters do not guarantee compatibility.

4. Do not force the plug into an outlet.
5. If the plug does not fit easily into the outlet, discontinue use.
6. Never wrap the cord tightly around the appliance. Never allow the power supply cord to be pulled, twisted, or severely bent. Untwist the cord before use.

DANGERS-WARNINGS

To reduce the risk of burns, electrocution, fire or injury to persons:

1. Toothbrush can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Cleaning and user maintenance shall not be made by children without supervision.
2. This product is designed to operate within a range of 100-240 volts.
3. The supply cord cannot be replaced. If the cord is damaged the appliance should be scrapped.
4. The battery must be removed from the appliance before it is scrapped;
5. The appliance must be disconnected from the supply mains when removing the battery;
6. The battery is to be disposed of safely.
7. Use this product only for its intended use as described in this manual. This product is designed to clean your teeth, gums and tongue only. Do not use any attachments or brush heads other than those recommended by the manufacturer

IMPORTANT SAFEGUARDS

8. The toothbrush is a personal care appliance and is not intended for use on multiple patients in a dental office or institution.
9. If your toothpaste includes peroxide, baking soda or other bicarbonate (common in whitening toothpastes), thoroughly clean your brush head and handle to avoid cracking the plastic.
10. Never operate this product if it has a damaged cord, plug or charger.
11. Do not use this product if it is not working properly, if it has been dropped or damaged or dropped into water. Immediately discontinue use if the product appears damaged in anyway.
12. Keep the toothbrush, charger base and power cord away from heated surface.
13. If you experience discomfort or pain in your mouth, discontinue use of this product and contact a dentist.
14. Do not attempt to remove the battery unless the product life is end.
15. Do not clean the toothbrush in the dishwasher.
16. Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.
17. This appliance has a polarized plug (one blade is wider than the other). As a safety feature, this plug will fit in a polarized outlet only one way. If the plug does not fit fully in the outlet, reverse the plug. If it still does not fit, contact a qualified electrician. Do not attempt to defeat this safety feature.

18. Do not leave appliance when plugged in. Unplug from outlet when not in use and before servicing.
19. Do not pull or carry by cord, use cord as a handle, close a door on cord, or pull cord around sharp edges or corners. Do not run appliance over cord.
20. Do not handle plug or appliance with wet hands.
21. Turn off all controls before unplugging.

MEDICAL-WARNINGS

1. If you have had oral surgery or gum surgery within the past months, consult your dental provider before using the product.
2. Consult your dental provider if excessive bleeding occurs when using this product or minor bleeding continues to occur after one week of use.
3. Consult your physician or implant device manufacturer prior to using this product if you have questions concerning the combined use of this product with a pacemaker or other implanted medical device(s).
4. Consult your physician prior to using the toothbrush if you have any medical concerns.
5. Always unplug this product immediately after using.
6. Never use while sleeping or drowsy.

These instructions shall also be available in an alternative format, e.g. on a website


SAVE THESE INSTRUCTIONS.

GETTING STARTED

Charging the toothbrush

Be sure to turn off the toothbrush before charging.

Connect the USB cable to power supply, the charging indicator will flash, indicating the product is charging. After the battery is fully charged, the charging indicator will turn off.

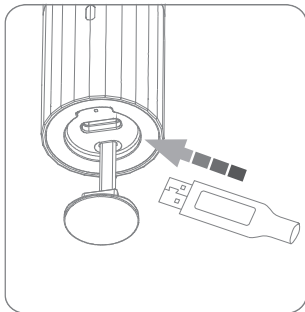

Notes:

1. It takes at least 4-6 hours to charge the battery fully. When fully charged, the toothbrush can be used 20 pre-timed brush cycles, or approximately 10 days, if brushing twice daily.
2. The toothbrush must be charged for 4-6 hours, if:
 - a. This is the first time using the toothbrush;
 - b. The toothbrush has not been used for more than three months;
3. To keep the battery fully charged, we recommend that you keep the toothbrush in the charger when not in use.
4. During charging, the toothbrush will get warm, but this is normal.
5. The operation duration will become shorter as the battery gets ageing. If the toothbrush operates only for one short time after being fully charged, this means that the battery is close to the end of its life cycle. Contact service center for support.

HOW TO USE

1. Slide the brush head onto the handle. Push brush head down until it locks into place on handle.

Note: There is small gap between color ring and the handle. (Fig 1)

2. Wet the bristles and apply a small amount of toothpaste.
3. Place the bristles against the teeth at a slight angle towards the gum line. (Fig 2)
4. Press the mode button to select mode. If you would like to keep the mode used last time, you need not select mode again.

a. CLEAN: Standard mode for superior teeth cleaning.

b. SOFT: Gently caress, carefully clean sensitive mouth.

c. WHITEN: Teeth cleaning with polishing and whitening.

d. MASSAGE: Pulsates for gentle gum stimulation.

e. POLISH: Polishes teeth to bring brightness naturally.

Note: During brushing, if pressing Mode button within 5s, toothbrush will work at next mode and start another 2 minutes cycle, if pressing after 5s, toothbrush will stop and the mode indicator will turn off.

5. Apply light pressure and gently move the brush head slowly across the teeth in a small back and forth motion.
6. In order to maximize cleaning effectiveness, divide your mouth in 4 sections. Start with Section 1 and brush each section for 30 seconds. The brushing action will pause between each 30 second interval, indicating when to move on to the next section. (Fig 3)
7. The toothbrush will automatically stop after 2 minutes of use.
8. Pull brush head out from the toothbrush. (Fig 4)
9. Rinse both the brush head and handle under running water and wipe them dry. (Fig 5)
10. Put the handle back into the charger. If you brush your teeth hard, the pressure sensor lamp will light up.

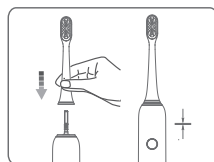


Fig 1

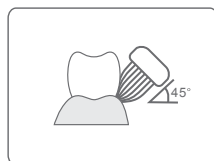


Fig 2

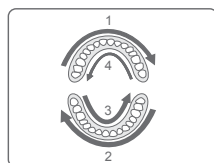


Fig 3

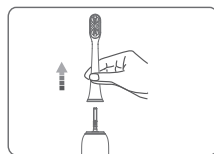


Fig 4

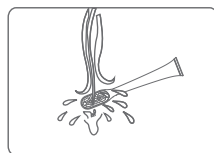


Fig 5


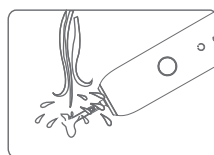



Fig 6

HOW TO USE

Take down the brush head from handle

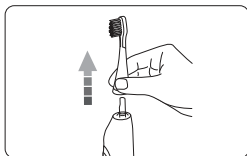

Rinse both the brush head and handle
under running water wipe them dry

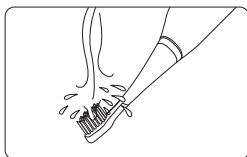


Put the handle back into the charger

MAINTENANCE

Be sure to remove the plug from the electrical outlet before cleaning the product.

Toothbrush


1. Remove the brush head.
2. Wipe all surface areas with a damp cloth and place the handle back into the charger.
3. Do not submerge the handle into water. Rinse the brush head and bristles after each use.
4. Do not clean in the dishwasher.

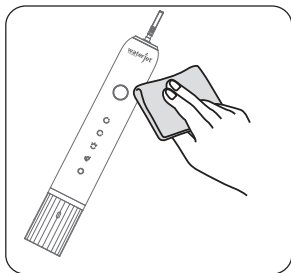

Charger

1. Using a soft cloth, wipe off any dirt on the charger and the storage.
2. Clean the plug periodically to ensure it remains in good working order.

The accumulation of dust or dirt on the plug can cause an electric shock, overheat or fire.

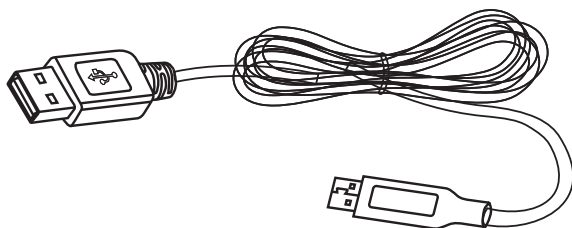


DIPOSAL

The power source of this product is Lithium-ion battery which can be recycled. At the end of the product's life, always remove the battery before you discard the product, and recycle or dispose of the battery according to local waste management requirements. If you have recycling questions, contact your local waste management office.


Removing and recycling the battery

Do not disassemble the product at any time other than when disposing of it. Doing so earlier will destroy the product. Observe basic safety precautions. The appliance must be disconnected from the supply mains when removing the battery.


1. Before removing the battery, be sure to completely drain the unit of power.

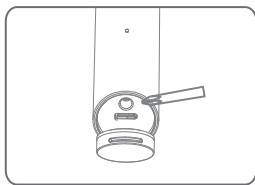


2. Use a screwdriver to remove the screw and take the internal components of the handle out.

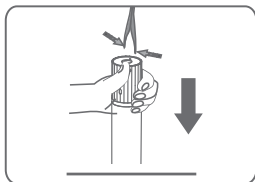

3. Remove the screw that fixes the bracket and PCBA, then melt the tin and take the battery out.

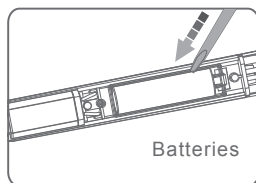

4. Recycle or dispose of the battery properly. Ask your local waste management office if they accept the remaining toothbrush components for recycling.

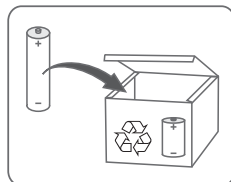

Be careful not to short-circuit the negative and positive poles of the battery. Do not connect the appliance to the mains again after you have removed the battery.

TROUBLESHOOTING

PROBLEM	POSSIBLE CAUSE	SOLUTION
The toothbrush doesn't work.	You have just purchased the product or it has not been used for more than 3 months.	Charge the main unit for at least 4 hours.
Weak or diminished vibration.	The toothbrush is pressed hard against your teeth.	Place the brush head gently against your teeth.
The toothbrush operates only for a short period of time even immediately after being charged.	The battery has reached the end of its life cycle.	Contact service center for support.
	Charging time was insufficient.	Charge the main unit for at least 4 hours.

WARRANTY

WaterJet will repair or replace (at our option) your unit free of charge for 12 months from the date of purchase if the appliance is defective in workmanship or materials.

This guarantee does not cover: damage due to improper use, normal wear or use as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorized persons and if original WaterJet parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your purchase receipt to an authorized WaterJet distributor/customer service centre in your country. Purchaser is responsible for all postage and freight charges and any other expenses incidental to claiming against this Warranty.



WaterJet

600 W Santa Ana Blvd Ste 114A PMB 224

Santa Ana, CA 92701, USA

Made in China

www.waterjet-flosser.com